

# SIZEWISE™: Alternation Therapy



• **SW Alternate™**

• **Mighty Air™**  
Bariatric

• **Consistency of design & operation facilitates staff in-services on product.**


**REAL THERAPY for REAL PEOPLE™**

**The Sizewise alternation therapy mattress systems** are therapeutic air support surfaces that are used for the prevention and treatment of pressure ulcers. Sizewise alternation therapy offers active pressure redistribution in two therapy modes: 1 in 3 Alternation Low Air Loss and Static Low Air Loss. Alternation low air loss therapy is for the treatment or prevention of skin breakdown, such as stages 1-4 pressure ulcers as well as for patients with limited mobility and respiratory issues. The 1 in 3 alternation system means every 3rd cell deflates depending on the programmed cycle time (between 3 to 20 minutes).

Static low air loss keeps skin cool and dry with continuous “true low air loss” which means constant cool air blows through the vapor permeable top covers that are especially designed to eliminate friction and shear and draw moisture away from the patient. This is ideal for your moisture management protocols. Under this modality, the caregiver can adjust the softness and firmness of the mattress based on the patient’s comfort or clinical decision.

The SW Alternate™ is used as a replacement mattress and can be placed on any standard size hospital bed frame including Hill-Rom® step-deck® frames. It can support patients weighing up to 600 pounds. The SW Mighty Air™ (bariatric) is also used as a replacement mattress in conjunction with the Bari-Rehab Platform™ and supports patients weighing up to 1,000 pounds.

## 1 in 3 Alternation Therapy

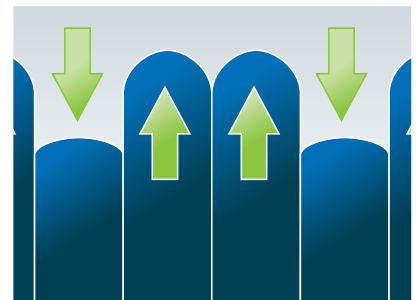

1740 liters of air per minute for maximum microclimate control



800.814.9389 | [sizewise.net](http://sizewise.net)

# SIZEWISE™: Alternation Therapy

## SW ALTERNATE™ SPECS

### CONTROL UNIT

Liters of Air Per Minute.....	1740
Alternating Cycle Time.....	Programmable from 3 to 20 minutes
Zoned Low Air Loss.....	Yes
Pulsating Low Air Loss.....	N/A
Continuous Lateral Rotation.....	N/A
Alarms.....	Power Failure & Low Pressure
Warranty.....	1 year
Pressure Sensor Technology.....	Yes
Patented CPR Release.....	Yes
Dimensions.....	18" W x 7 ½" D x 11 ½" H
Unit Weight.....	16 lbs.

### MATTRESS

Weight Capacity.....	600 lbs. in Static
Mattress Size.....	Standard 35"
Warranty.....	2 Years Limited*
Top Cover Material.....	70 Denier Polyurethane Coated Nylon Taffeta
Height Inflated.....	Approx. 10"
Mattress Weight.....	37 lbs.
Top Cover Warranty.....	90 days

### ELECTRICAL

Rated Voltage.....	110 volts
Rated Frequency.....	60 Hz
Rated Input Power.....	120 volts
Fuse Rating.....	125 volts 5 amps
Power Failure Alarm.....	Yes

## MIGHTY AIR™ SPECS

### CONTROL UNIT

Liters of Air Per Minute.....	1740
Alternating Cycle Time.....	Programmable from 3 to 20 minutes
Zoned Low Air Loss.....	Yes
Pulsating Low Air Loss.....	N/A
Continuous Lateral Rotation.....	N/A
Alarms.....	Power Failure & Low Pressure
Warranty.....	1 year
Pressure Sensor Technology.....	Yes
Patented CPR Release.....	Yes
Dimensions.....	18" W x 7 ½" D x 11 ½" H
Unit Weight.....	16 lbs.

### MATTRESS

Weight Capacity.....	1,000 lbs.
Mattress Size.....	39" & 48" widths
Warranty.....	2 Years Limited*
Top Cover Material.....	70 Denier Polyurethane Coated Nylon Taffeta
Height Inflated.....	Approx. 10"
Mattress Weight.....	37 lbs.
Top Cover Warranty.....	90 days

### ELECTRICAL


Rated Voltage.....	110 volts
Rated Frequency.....	60 Hz
Rated Input Power.....	120 volts
Fuse Rating.....	125 volts 5 amps
Power Failure Alarm.....	Yes

\* Air Support limited warranty includes two (2) years on air cell welds; one (1) year on bottom cover, manifold and foam base; Ninety (90) Days – Top cover. All parts found defective within that period shall be repaired and/or replaced.

**MADE IN THE USA**


## other available THERAPIES

### Pulsation Therapy – SW Pulsate™

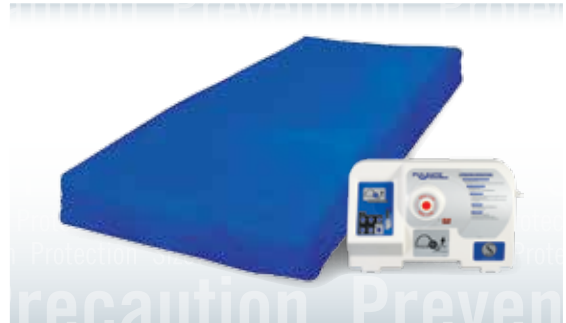

Bariatric version also available

### Rotation Therapy – SW Rotate™

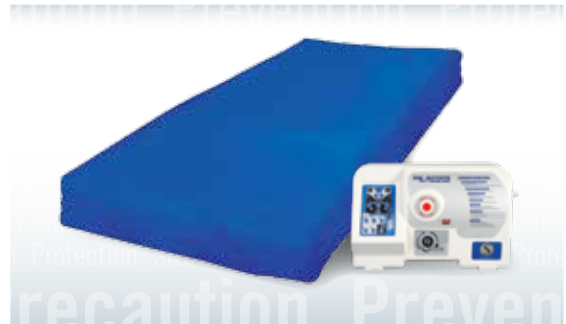

### Bariatric Rotation Therapy – Big Turn™



**SIZEWISE**

800.814.9389 | [sizewise.net](http://sizewise.net)