


SIZEWISE™: OR Table Pad™

A part of the **Sizewise Envy™ Line** – the industry's first surfaces that reduce exposure to harmful chemicals of concern.

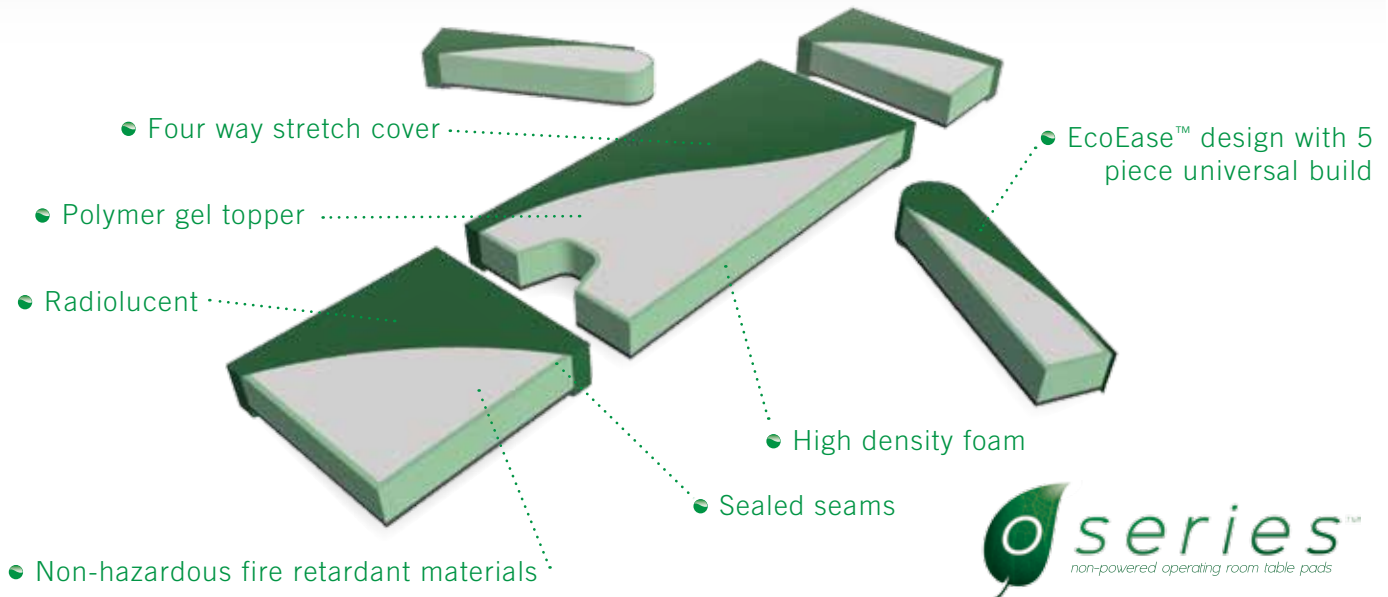

Sizewise Envy Line™ O Series™ is the industry's first non-powered, self-adjusting gel operating room table pad featuring F³ Free™ construction that reduces chemicals of concern in the fabric, foam & fire barrier. The five piece universal design is radiolucent and fits most OR tables. The top cover is a four way stretch cover that allows immersion into the pad for optimal pressure relief. Wipe down sterilization is made possible by the sealed seam cover assembly which is free from sewn stitch holes. The cover is fluid repellent and stain-resistant. Specially formulated polymer gel topper contours to the patient's body while providing low interface pressure to the patient. A high density foam base provides additional stability while enhancing comfort and pressure relief. EcoEase™ designed pieces are easily replaceable by section creating an economical and environmentally conscious product. This product meets the Practice Greenhealth Environmental Specifications for Healthcare Support Surfaces. Weight capacity is 500 lbs.

Multiple Sizes Available

Sizewise O Series™

Head Section	10-1/2" x 20" x 4"
Torso Section	42-3/8" x 20" x 4"
Leg Section	18-1/2" x 20" x 4"
Arm Board Pad	26" x 6" x 4"
Arm Board Riser	26" x 6" x 2"

Sizewise O Series™

Head Section	9" x 20" x 4"
Torso Section	40-1/4" x 20" x 4"
Leg Section	23-1/4" x 20" x 4"
Arm Board Pad	26" x 6" x 4"
Arm Board Riser	26" x 6" x 2"



800.814.9389 | sizewise.net

REAL THERAPY for REAL PEOPLE™

SW.11.04.11



created with recycled material & soy based ink

All specifications and equipment are subject to change without notice. Sizewise reserves the right to make improvements from time to time. Photos and drawings are representative of the products and may vary slightly from actual production models. Some items shown may include optional equipment.