

SIZEWISE™: Rotation Therapy



• **SW Rotate™**

• **Big Turn™**
Bariatric


• **Consistency of design & operation facilitates staff in-services on product.**

REAL THERAPY for REAL PEOPLE™

The Sizewise rotation therapy mattress systems provide continuous low air loss therapy in all function settings. If there are times when rotation is not desired, the pulsation mode can be used to provide active pressure redistribution by increasing and decreasing pressure over the entire body at preset cycle times. In addition, the optional respiratory therapy module allows mechanical percussion and vibration for the immobile patient. The control settings allow for three areas of the lung to be percussed and vibrated independently with varying degrees of intensity, from very soft percussion and vibration to a very aggressive degree of therapy. The nurse can begin the therapy and walk away knowing the control will automatically stop the therapy at a prescribed time.

The SW Rotate™ is used as a replacement mattress and can be placed on any standard size hospital bed frame including Hill-Rom® step-deck frames. It can support patients weighing up to 600 pounds. The Big Turn™ (bariatric) used in conjunction with the Bari-Rehab Platform™ supports patients weighing up to 1,000 pounds.

Turn Assist Feature

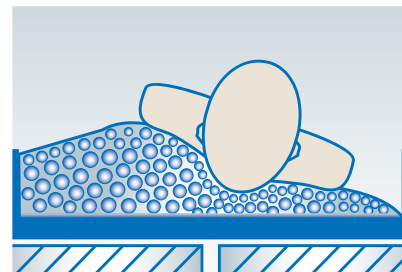

1740 liters of air per minute for maximum microclimate control



800.814.9389 | sizewise.net

SIZEWISE™: Rotation Therapy

SW ROTATE™ SPECS

CONTROL UNIT

Liters of Air Per Minute.....	1740
Alternating Cycle Time.....	N/A
Zoned Low Air Loss.....	Yes
Pulsating Low Air Loss.....	3, 5, or 10 minutes
Continuous Lateral Rotation.....	Programmable 1 - 30 minutes
Rotation Training Mode.....	Yes
Alarms.....	Power Failure & Low Pressure
Warranty.....	1 year
Pressure Sensor Technology.....	Yes
Patented CPR Release.....	Yes
Dimensions.....	18" W x 7 ½" D x 11 ½" H
Unit Weight.....	16 lbs.

MATTRESS

Weight Capacity.....	350 lbs. in Rotation/ 600 lbs. in Static
Mattress Size.....	Standard 35"
Percussion Modules.....	Add-on Feature
Warranty.....	2 year
Height Inflated.....	Approx. 10"
Mattress Weight.....	37 lbs.
Top Cover Material.....	70 Denier Polyurethane Coated Nylon Taffeta
Top Cover Warranty.....	90 days

ELECTRICAL

Rated Voltage.....	110 volts
Rated Frequency.....	60 Hz
Rated Input Power.....	120 volts
Fuse Rating.....	125 volts 5 amps
Power Failure Alarm.....	Yes

BIG TURN™ SPECS

CONTROL UNIT

Liters of Air Per Minute.....	1740
Alternating Cycle Time.....	N/A
Zoned Low Air Loss.....	Yes
Pulsating Low Air Loss.....	3, 5, or 10 minutes
Continuous Lateral Rotation.....	Programmable 1 - 30 minutes
Rotation Training Mode.....	Yes
Alarms.....	Power Failure & Low Pressure
Warranty.....	1 year
Pressure Sensor Technology.....	Yes
Patented CPR Release.....	Yes
Dimensions.....	18" W x 7 ½" D x 11 ½" H
Unit Weight.....	16 lbs.

MATTRESS


Weight Capacity.....	1,000 lbs.
Mattress Size.....	39", 48"
Percussion Modules.....	Add-on Feature
Warranty.....	2 year
Height Inflated.....	Approx. 10"
Mattress Weight.....	37 lbs.
Top Cover Material.....	70 Denier Polyurethane Coated Nylon Taffeta
Top Cover Warranty.....	90 days

ELECTRICAL

Rated Voltage.....	110 volts
Rated Frequency.....	60 Hz
Rated Input Power.....	120 volts
Fuse Rating.....	125 volts 5 amps
Power Failure Alarm.....	Yes


other available THERAPIES

Pulsation Therapy – SW Pulsate™

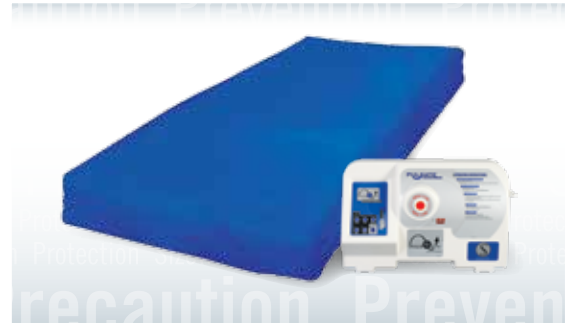

Bariatric version also available

Alternation Therapy – SW Alternate™

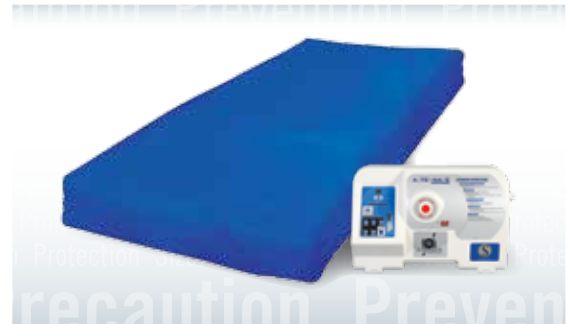

Bariatric Alternation Therapy – Mighty Air™



800.814.9389 | sizewise.net