

SIZEWISE™: Evolution™ Expandable

Capital Purchase Only

• Ergonomic Design

• Expandable deck 39", 42", 48" width

• Lock & Steer System

• Independent Lockout Features

• Rest Secure System™

REAL THERAPY for REAL PEOPLE™

Introducing Evolution™ Expandable – the next generation in versatile, patient care beds. With an 850 pound capacity and class leading 9" low position, the Evolution™ Expandable is a safe option for 'at risk fall' patients. The deck height of 30.75" allows for an adjustable ergonomic working environment for care providers of all sizes and the bed can transport a patient in any high or low position. The deck width expands from 39" to 42" to 48" allowing for flexibility in meeting patient needs. When combined with our therapeutic support surfaces, the Evolution™ Expandable can provide alternation therapy, pulsation therapy or rotational therapy with percussion and vibration to accommodate a wide array of patient needs. The Rest Secure System™ is an integrated scale, bed exit alarm and protocol timer that interfaces with the nurse call system. Caregiver lockout controls, nurse call interface, bed exit alarm, cardio chair positioning, 20 degree Trendelenberg and Reverse Trendelenberg make the Evolution™ Expandable the ideal bed for any care environment.

Cardio Chair



Trapeze Option



Rest Secure System™



Expandable Widths: 39", 42", 48"



SIZEWISE

800.814.9389 | sizewise.net

SIZEWISE™: Evolution™ Expandable

available THERAPIES

BED

Low Position - Deck to Floor	9"
High Position - Deck to Floor	30.75"
Sleep deck width	Expandable 39", 42", 48"
Sleep deck length	82"
Maximum overall width w/ rails	52"
Maximum overall length w/o bumper	89"
Patient weight capacity	850 lbs

SIDE RAILS

Type	Fold Down
Included with purchase	One Left, One Right, Head Rail

HEAD/FOOT BOARDS

Included in Purchase	Yes
Type	High Impact Plastic

BED ANGLES

Head Section	0° - 65°
Knee Section	0° - 50°
Max Trend	0° - 20°
Max Reverse Trend	0° - 20°

CASTER & BRAKING SYSTEM

Caster Diameter	4"
Dual Brake/Steer System	4"

SCALE SYSTEM

Included in Purchase	Yes
Weigh in any position	Yes

ELECTRICAL SPECS

Emergency Battery back-up	Included
Safety Listing	EN60601
Foot End Bed Controls	Included
Power Requirements	120v, 60 cycle

MATTRESS OPTIONS

Low Air Loss Surface	Optional
Alternation Therapy	Optional
Full Body Lateral Rotation Therapy	Optional
Pulsation Therapy	Optional
Non-Powered Surface	Optional
Tri-Laminate Surface	Optional

OPTIONS

Trapeze	Yes
---------------	-----

PATENTS

U.S. Patent No. 7,013,510
U.S. Patent No. 7,334,277



Alternation Therapy – Sapphire 1100EC™




Pulsation Therapy – Sapphire Pulse™



Rotation Therapy – Comfort Turn™



Control Panel with Independent Lockout Features

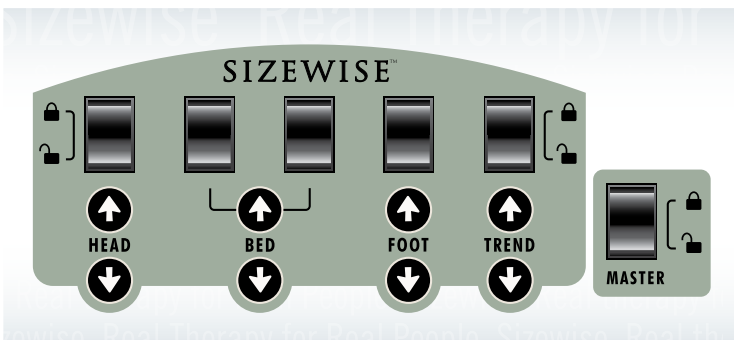

800.814.9389 | sizewise.net