



# SIZEWISE™: Stretcher Pad

A part of the **Sizewise Envy™ Line** – the industry's first surfaces that reduce exposure to harmful chemicals of concern.

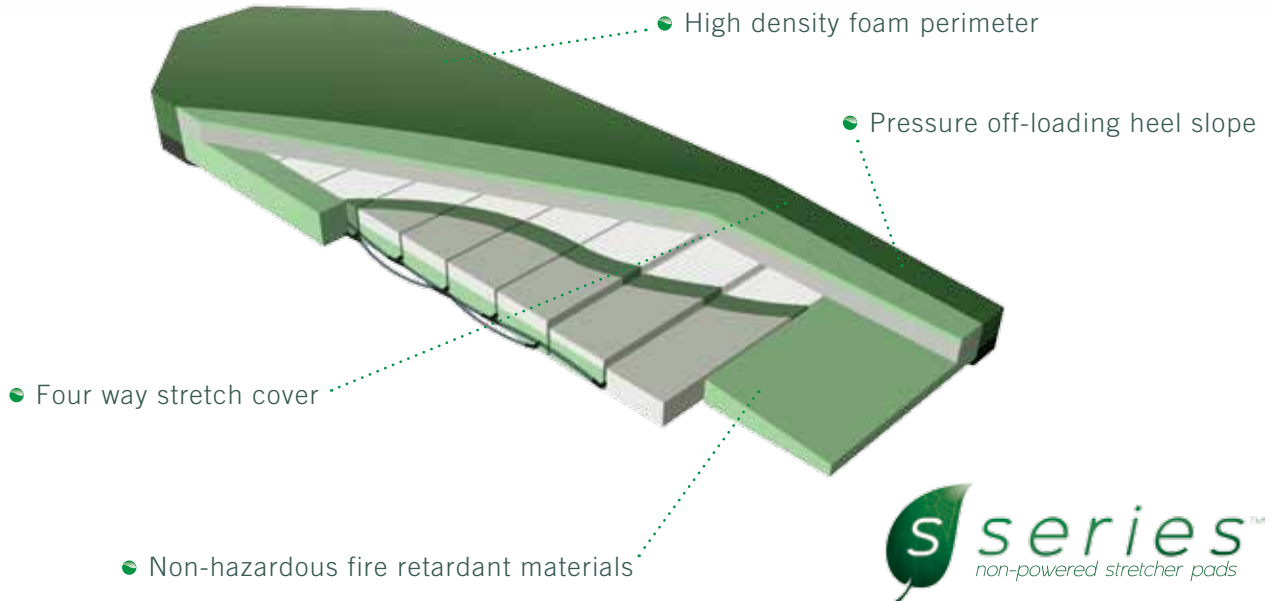

**Sizewise Envy Line™ S Series™** is the industry's first non-powered, self-adjusting gel stretcher pad system featuring F3 Free™ construction that reduces chemicals of concern in the fabric, foam & fire barrier. The top cover is a stain-resistant and fluid repellent four way stretch cover that allows the patient to immerse into the pad for optimal pressure relief. Specially formulated topper contours to the patient's body while providing low interface pressure to the patient. The heel slope is designed to off-load pressure from the heels while the perimeter of high density foam aids in transfer and provides additional stability while enhancing comfort and pressure relief. Velcro runs from the top of the pad to the bottom to reduce product migration when adhered to a stretcher. Custom construction is available in the S Series, including a model that consists solely of tri-laminate pressure redistribution foam. These products exceed the Practice Greenhealth Environmental Specifications for Healthcare Support Surfaces. Weight capacity is 500 lbs., bariatric options are available.


## Pressure Mapping - 50 mmHg Scale Shown

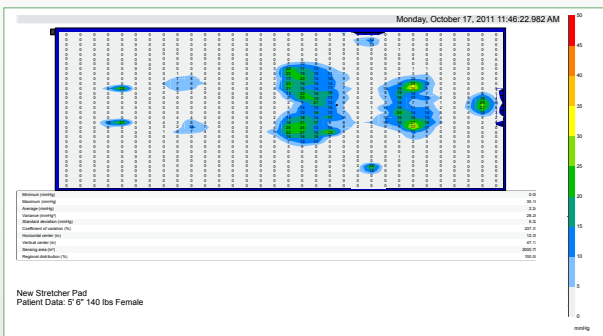

## S Series™ Non-Powered Stretcher Pad Featured

**MADE IN THE USA**

## Optional Model - Tri-Laminate Pressure Redistribution Foam

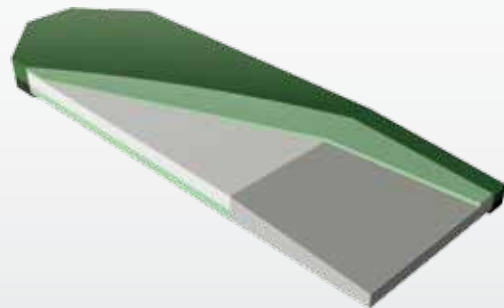

**S** SIZEWISE

800.814.9389 | [sizewise.net](http://sizewise.net)

REAL THERAPY for REAL PEOPLE™

SW.11.30.11



created with  
recycled material  
& soy based ink

All specifications and equipment are subject to change without notice. Sizewise reserves the right to make improvements from time to time. Photos and drawings are representative of the products and may vary slightly from actual production models. Some items shown may include optional equipment.