

SIZEWISE™: Bari-Rehab Platform²™



REAL THERAPY for REAL PEOPLE™

Introducing Bari-Rehab Platform²™ – setting a new state-of-the-art standard for bariatric care. With a 1,000 pound capacity and expandable width of 39" or 48" the Bari-Rehab Platform² is the ideal bed for all levels of care for the bariatric patient. This innovative bed design allows patients to become actively involved in their rehabilitation plan of care. Features such as the open side rail, low working height of the deck and simplified hand control all enhance the patient's independence as well as decrease the patient's reliance on the caregiver. When combined with our low air loss therapeutic support surfaces, the Bari-Rehab Platform² can provide alternation therapy, pulsation therapy, or rotational therapy with percussion and vibration to accommodate a wide array of patient conditions and caregiver needs. The integrated scale system delivers easy, accurate weights at the touch of a button. All models are fully electric and come standard with cardio chair positioning, up to 15° Trendelenburg and Reverse Trendelenburg, and include a built in Head Angle Indicator.

Head Angle Indicator



Power Drive Joystick



800.814.9389 | sizewise.net

SIZEWISE™: Bari-Rehab Platform²™

BED – Patent Number: 5,669,0190

Low Position - Deck to Floor	15.5"
High Position - Deck to Floor	30"
Expandable width.....	Yes
Sleep deck width.....	39" to 48"
Sleep deck length	80" or 86"
Maximum overall width w/ rails.....	42" or 51"
Maximum overall length w/o bumper	90" or 96"
Patient weight capacity.....	1000 lbs

MATTRESS – Patent Number: 7,225,489

Included in Purchase	Yes
Mattress Height.....	7"
Seamless Expandable Mattress	Included
Fire Safety.....	Cal 129
Fire Safety TB 603	Optional

SIDE RAILS

Type	Removable, 3 Position Fold Down
Material	Kydex®
Included with purchase.....	Two Sets

HEAD/FOOT BOARDS

Included in Purchase	Yes
Material	Kydex®
Accessory Bar.....	Yes
Info Pouch.....	Yes
Ergonomic	Yes

POWER DRIVE

Accepts Power Drive Option	Yes
----------------------------------	-----

BED ANGLES

Head Section.....	0° - 60°
Knee Section.....	0° - 50°
Max Trend.....	0° - 15°
Max Reverse Trend.....	0° - 15°

CASTER & BRAKING SYSTEM

Dual Brake/Steer System.....	5"
------------------------------	----

SCALE SYSTEM

Included in Purchase	Yes
Weigh in any position	Yes

(except Trend/Rev Trend)

ELECTRICAL SPECS

Battery back-up.....	Optional
Safety Listing.....	UL 544
Power Requirements.....	120v, 60 cycle

MATTRESS OPTIONS

Low Air Loss Surface.....	Optional
Alternation Therapy.....	Optional
Full Body Lateral Rotation Therapy	Optional
Pulsation Therapy	Optional
Non-Powered Surface	Optional
Tri-Laminate Surface.....	Optional

OPTIONS

Trapeze.....	Yes
--------------	-----



available THERAPIES


Alternation Therapy – Mighty Air™



Rotation Therapy – Big Turn™



Pulsation Therapy – SW Pulsate™

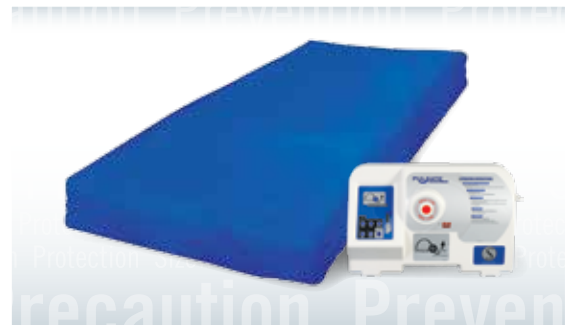

800.814.9389 | sizewise.net